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breo

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EN

Things Prohibited

	Read the Instructions Manual carefully. Avoid long-time continuous use of this product to prevent muscle pain caused by excessive massage.
	Do not use this product while driving or conducting precision operations.
	Do not use on surgical wounds or damaged/broken skin or area with inflammation or congestion.
	Do not use if suffering from heart disease (or with a history of heart disease), acute disease, skin disease, malignant tumors, blood disease (such as bleeding), osteoporosis, soft tissue disease, stroke or other physical discomforts.
	Not for use by persons who feel painful to move or with physical or mental discomforts or intellectual disabilities. In case of using, follow the doctor's advice.

3

	Do not use if pregnant or in menstrual period. In case of using, follow the doctor's advice.
	Not for use by children under 8. Children aged from 8 to 12 are not recommended to use this product unless they have been given supervision or instruction concerning use of this product by parents or guardians.

Caution

	This is a household electronic healthcare product and shall not be used as medical device or for diagnosis, treatment or other commercial purposes.
	Before use, remove accessories around the massage area or other items that may affect the massage.
	Stop using if pain or discomfort occurs in the massage area.

4

	During use, if exceptions occur to the product, cut off power supply immediately and contact the customer service hotline. Do not disassemble or repair it by yourself.
	Persons who are not sensitive to heat shall be careful during use as the surface of the product may be heated.
	When using under other specific conditions, please observe safety regulations.

Product Care

	Clean the surface of the product with soft wet cloth or a 75% alcohol pad. Do not rinse.
	Charge the product with the charging device provided by our Company. Beware of electricity safety.
	Keep the product away from water, fire and corrosive environment.

EN

NOTES

Must Know

This massage head is made of food-grade silicone/metal that can directly contact the skin. Please read the Instructions Manual carefully before use to avoid skin damage.

※Usage Frequency

Twice daily (morning and evening) for no longer than 10 minutes each is recommended.

※ Cautions before Use

- 1 Do not use with wet hands.
- 2 Stop use and seek medical attention if pain or burning occurs in the massage area.

※ After Use

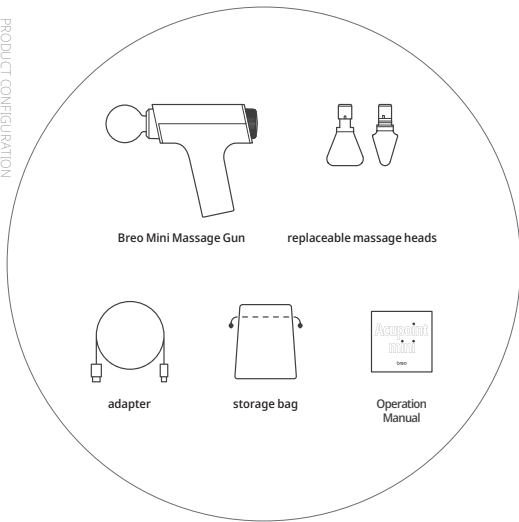
The massage area may become red and hot (warm) in a short time, which is normal and will recover soon.

※ How to Clean

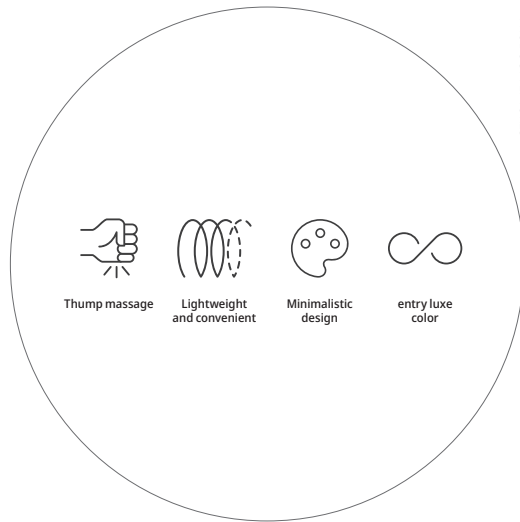
- 1 Clean the surface of the product with soft wet cloth or a 75% alcohol pad. Do not rinse. (For plastic and other smooth surface)

※ Others

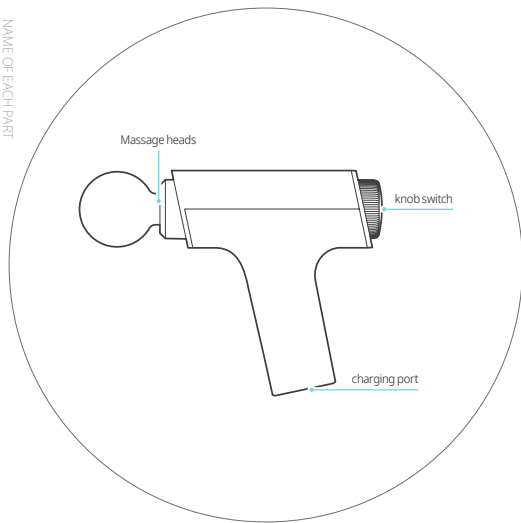
This product may drone, purr or rap during normal use, which is normal.

7

8

EN

Instructions

Indicator

Red light blinks when battery is low (slow blinking: once every 1s).

Red light blinks when charging.

Green light on once fully charged.

Red light on in case of overcurrent protection, and the device stops working.

The indicator turns purple when the pressure returns to normal, and the device continues working.

Switch 0

Turn off the device.

Switch 01

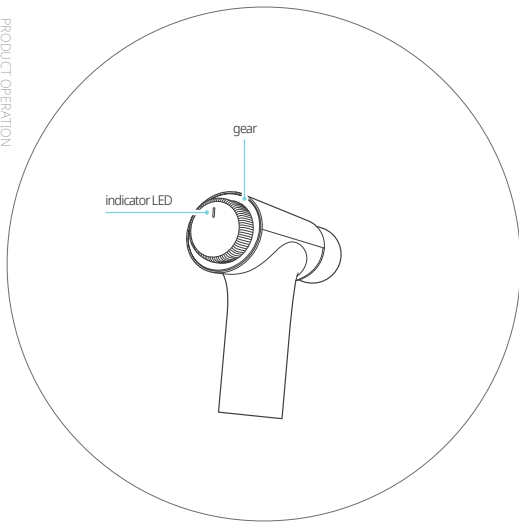
Gentle massage mode.

Switch 02

Moderate massage mode.

Switch 03

Strong massage mode.moving and rotating from back to front with heat function off.

11

EN**Charging Methods**

1 Insert the DC plug into the DC port of the device.

2 Plug the power adapter into the socket.

Red light blinks when charging and keeps on once fully charged.

※It is normal for the device and the power adapter to heat up a little while charging.

Charging Time

1 The charging time is 2 to 3 hours.

2 Once fully charged, it can be used 15 times for 10 minutes each time.

(It may vary in environments of different temperatures)

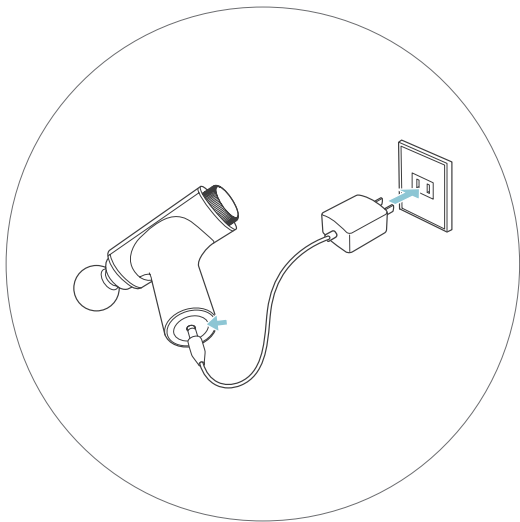
Battery

The indicator blinks when the battery is low. Please charge before use again.

After Charging

Unplug the power adapter and the DC plug.

12

13

EN

How to Use

Do not use while charging. Please unplug the power adapter once fully charged.

Power on

rotate the knob switch to turn on the power, and the device starts with the indicator on.

Mode description

Adjust your favorite massage mode by rotating the knob switch according to your own needs.

Switch 0

Turn off the device.

Switch 01

Gentle massage mode, body refreshing (2,000 rpm)

Switch 02

Moderate massage mode, relax the body (2,600 rpm)

Switch 03

Strong massage mode, relieving painful and sore muscles (3,100 rpm)

14

Power off

Automatic shutdown: the device will automatically shutdown with the indicator off after using for a total of 10 minutes. Turn the knob switch to OFF manually and restart if needed.

Manual shutdown: turn the knob switch to O, and the device will turn off with the indicator off.

Protection mechanism

1 Equipped with overcurrent protection. Excessive pressure on the massage head may affect health. To ensure safety, the device will automatically stop working, with the red indicator on. In this occasion, please move away the massage head from the body until the indicator turns purple before returning to normal use.

2 Equipped with over-heat protection. Continuous use may cause the motor in the handle to heat up. The device will automatically stop working when it detects an over-heat in the motor, with the red light on. In this occasion, please stop using the device until the indicator turns purple.

EN

Operation steps for battery removal

- 1 Cut off the external power supply, and remove the product housing;
- 2 Disconnect the power supply to the battery;
- 3 Remove the battery;
- 4 Place the waste battery at a fixed recycling station;
- 5 Please observe local environmental protection laws and waste disposal guidance and protect the environment.

Notes

- Remove the battery out of the appliance before scrapping the appliance;
- The appliance must be disconnected from power supply before removing the battery;
- The battery should be disposed of safely.

Except professionals, no one shall attempt to remove the battery!

EN Troubleshooting

If any faults occur during use, please check the following. If problem still exists or the product doesn't work well, please cut off power supply and send it to professional and technical staff for repair. Do not disassemble it by yourself!

Symptom	Possible cause	Basic solution
Power-on failure	Check whether the battery level is too low	If yes, please charge the battery timely
Low massage intensity	Check whether the battery level is low	If yes, please charge the battery timely

If you have any issues, questions or feedback, please send an email to support@breocare.com or call us at 866-546-3587. We are always happy to help.

EN

PROTECTION

Storage

When the device is stored for more than 3 months, check the device and recharge it before using. Do not store the device closely against the ground, wall or ceiling.

Store in a well-ventilated place and avoid direct sunlight and harmful erosion.

Transport

The product is packed with simple anti-vibration measures for transportation by air, railway, highway and sea. Avoid rain, snow and collision shall be avoided and do not leave it upside down.

Power Check

Check the power supply. Do not start the machine if power supply exceeds the allowed input.

Maintenance

Take measures against moisture, rodents and pollution. Turn the machine off before cleaning. Store the machine in a safe place after use. Do not place heavy items on the machine. If any fault occurs, call professional and technical staff for repair. Do

not disassemble the machine if you are non-professional.

Disposal and Recycling

This product cannot be discarded as household waste. Please follow local environment laws and guidelines when you have decided to discard this product.

Please remember to respect the local regulations: hand in the non-working electrical equipments to an appropriate waste disposal centre. The packaging material is recyclable. Dispose of the packaging in an environmentally friendly manner and make it available for the recyclable material collection-service.

This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

This device complies with Part 15 of the FCC Rules. Op-

eration is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. The statements should be displayed in the user manual: changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

EN

Product name	Acupoint Mini
Product model	Sport-M0210
Rated voltage	15V ---
Rated power	25W
Battery capacity	1900mAh
Product weight	455g
Product size	L142.5×W94.8×H48(mm)

Shenzhen Breo Technology Co., LTD

19/F, Hisense Southern Building, No.1777 Chuangye Road, Yuehai Street,
Nanshan District, Shenzhen 518054, China

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成品尺寸: L85 x W85mm 展开尺寸: L85 x W170mm ©copyright:zhao 20200810 材质: 白色道林纸